# Detection of pyrogens in fetal bovine serum (FBS) with the PyroMAT<sup>™</sup> System

### Introduction

### What is a pyrogen?

A pyrogen is, by definition, a substance that produces a rise in temperature in a human or animal. Pyrogens constitute a heterogeneous group of contaminants comprising microbial and non-microbial substances. The most widely known pyrogen is the endotoxin (LPS = Lipo-Polysaccharide), which is produced by gram-negative bacteria. Other microbial substances include those derived from gram-positive bacteria like Lipoteichoic Acid (LTA), particles from viruses and pyrogens originating from yeasts and fungi. Nonmicrobial pyrogenic substances can be rubber particles, microscopic plastic particles or metal compounds in elastomers.

### Why to conduct a pyrogen test?

Pyrogenic substances in pharmaceutical products can induce life-threatening fever reactions after injection into the human body. Therefore, it is a regulatory requirement to test such products for pyrogens to ensure product quality and patient safety.

Purpose of the test is to prove that the amount of pyrogens contained in the product will not exceed a certain threshold, known as the contaminant limit concentration (CLC), that will guarantee the patient safety.

The monocyte activation test (MAT) method has been qualified and validated for the detection of pyrogens by the European Center for the Validation of Alternative Methods (ECVAM) in 2005 and by the Interagency Coordinating Committee on the Validation of Alternative Methods (ICCVAM) in 2008.

It has been among the compendial methods for pyrogen detection in the European Pharmacopeia since 2010 (Chapter 2.6.30) [1].

The MAT is also mentioned by the FDA "Guidance For Industry – Pyrogen and Endotoxins testing: Questions and Answers" as an alternative to the rabbit pyrogen test which should be validated according to USP <1225>. Additionally, the USP <151> Pyrogen Test mentions that, "A validated, equivalent *in vitro* pyrogen or bacterial endotoxin test may be used in place of the *in vivo* rabbit pyrogen test, where appropriate."

### **Principle of the MAT**

The monocyte activation test (MAT) is the human *in vitro* alternative to the rabbit pyrogen test, and allows the detection of the full range of pyrogens, including endotoxins and non-endotoxin pyrogens (NEPs).

By putting the product to be tested in contact with human monocytic cells, it will mimic what happens in the human body: in presence of pyrogens, the monocytes are activated and produce cytokines such as interleukin-6.

The cytokines are then detected using an immunological assay (ELISA) involving specific antibodies and an enzymatic color reaction.





### Principle of the PyroMAT<sup>™</sup> System

The PyroMAT<sup>™</sup> System uses cryo-preserved Mono-Mac-6 (MM6) human monocytic cells as a source of monocytes.

The response to pyrogenic substances is determined by measurement of Interleukin-6 (IL-6) produced by the Mono-Mac-6 cells. For this purpose, the ELISAmicroplate supplied in the kit is coated with an antibody specific to IL-6.

IL-6 molecules released by MM6 cells supernatant during incubation phase are transferred in the ELISA plate, and bound by the immobilized primary antibody.

A secondary antibody, linked to an enzyme, is added to form an IL-6 bound complex. After washing any unbound molecules, the IL-6 bound complex is detected in a color reaction started by the addition of an appropriate substrate.

The color development is proportional to the amount of initial IL-6 production in the supernatant and measured with an absorbance reader.

### Quantification of pyrogens with the MAT

For the quantification of the pyrogenic load of a sample, method A in accordance with the European Pharmacopeia [1] can be conducted.

Method A involves a comparison of the preparation being examined with a standard endotoxin dose-response curve. The contaminant concentration of the preparation being examined is to be less than the CLC (Contaminant Limit Concentration) to pass the test.

To ensure both the precision and validity of the test, preparatory tests need to be conducted to assure that:

- The criteria for the endotoxin standard curve are satisfied
- The solution does not interfere with the test
- The test detects endotoxin and non-endotoxin contaminants
- The solution does not interfere with the detection system

### **Material and Equipment**

To perform the MAT and a product specific validation, we recommend using:

- PyroMAT<sup>™</sup> Cells (Ref: Pyr0MATCELLS)
- PyroMAT<sup>™</sup> Kit (Ref: Pyr0MATKIT)
- Reference Standard Endotoxin (Ref: 1.44161.0001).
- NEP Control HKSA (Ref: MATHKSA)
- NEP Control Flagellin (Ref: MATFLAGELLIN)
- IL-6 control (Ref: Pyr0MATIL6)

Additional equipment and consumables required:

- Incubator, 37 °C, humidified
- Water bath (37 °C)
- Microplate reader to measure absorbance at 450 nm and 630 nm (reference wavelength)
- Cryo-freezer (< -80 °C)</li>
- Freezer (-20 °C)
- Refrigerator (2-8 °C)
- Vortexer
- 50 mL centrifuge
- Multichannel pipettes with suitable containers
- Adjustable pipettes: (10  $\mu L$  100  $\mu L$ ; 100  $\mu L$  1000  $\mu L$ ) with suitable sterile, pyrogen-free pipette tips
- Pyrogen-free glass tubes
- 2 mL endotoxin-free reaction tubes with caps
- 50 mL endotoxin-free tubes with caps



### Quantification of pyrogens with the PyroMAT<sup>™</sup> system

### Preparatory tests – Product Specific Validation (PSV)

#### European pharmacopeia, chapter 2.6.30

Before routine testing of a pharmaceutical product with MAT, a product specific validation (PSV) must be performed according to the method chosen for routine testing to ensure the validity of the criteria for the endotoxin standard curve, the detectability of endotoxin and non-endotoxin contaminants in the sample and that the sample does not interfere with the test or detection system.

Interferences with the test can be removed by diluting the product up to a certain limit, referred to as the maximum valid dilution (MVD).

The MVD is the maximum dilution factor at which it is still possible to detect the pyrogen limit (i.e., the CLC). It is directly linked to the limit of detection (LOD) of the system.

The more sensitive the system is, the more the product can be diluted to remove interferences.

The MVD of a test solution can be calculated using the following formula:

$$MVD = \frac{CLC \times C}{LOD}$$

CLC = Contaminant Limit Concentration (EU/mg or EU/mL)C = Concentration of the test solution (mg/mL or mL/mL) LOD = Limit of Detection (EU/mL).

The CLC is the acceptance criterion for the pass/fail decision, expressed in endotoxin equivalents per milligram or milliliter (EEU/mg or EEU/mL) or per unit of the biological activity of the product.

It is calculated by the following expression:

$$CLC = \frac{K}{M}$$

K = threshold pyrogenic dose per kilogram of body mass (EU/kg)

M = maximum recommended bolus dose of product per kilogram of body mass (mg/kg or mL/kg).

When the product is to be injected at frequent intervals or infused continuously, M is the maximum total dose administered in a single hour period.

When testing for interfering factors, dilutions of the preparation being examined with geometric steps not exceeding the MVD should be performed. The same dilutions spiked with endotoxin at a justified concentration (in case of method A, a concentration near the estimated middle of the endotoxin standard curve) should then be performed. These dilutions must be tested in parallel in the same experiment and together with an endotoxin standard curve, which shall be used to calculate the concentration of endotoxin-equivalents in each solution.

The mean recovery of the added endotoxin spike is then calculated for each dilution. The test is considered free of interference when recovery of the added endotoxin is within the range of between 50 and 200%.

Where practicable, interference testing should be performed on at least 3 different lots of the preparation being examined to investigate on possible batch-tobatch variation. If the interference cannot be removed by dilution or specific sample preparation of the product within the MVD range, Method C is preferred over Method A and B.

For validation of the detection of non-endotoxin contaminants, historical batches that have been found to be contaminated with non-endotoxin contaminants causing positive responses in the rabbit pyrogens test or adverse drug reaction in man can be used. Where those batches are not available, validation should be done including at least 2 non-endotoxin ligands for toll-like receptors, one of which is to be spiked into the preparation being examined.

Once the optimum dilution of the preparation being examined has been identified, this dilution needs to be tested for interference in the detection system. The agreement between a dilution series in presence and absence of the preparation being examined is to be within  $\pm$  20% of optical density.

### Sample specifications: FBS (fetal bovine serum)

FBS derived from clotted blood is the most widely used undefined supplement in eucaryotic, especially mammalian, cell culture. Though there are efforts to establish human or synthetic alternatives, fetal bovine serum still is a product that often is used in the production of vaccines and therapeutics.

Being a raw material used within the production process of injectables administered to man, FBS has to be tested for harmful substances like endotoxin or other pyrogenic substances like  $(1 \rightarrow 3)$ - $\beta$ -d-glucan [2,3]. Quality of FBS is also important for researchers using cell cultures as pyrogenic contaminants may create problems by affecting the bioactivity of the cultured cells and therefore influence the experimental results obtained [4].

For the examined FBS the criteria of quality was an endotoxin level of equal or smaller than 10 EU/mL. Therefore, the contaminant limit concentration of the product is:

CLC = 10 EU/mL

For the PyroMAT<sup>™</sup> system, the LOD is 0.05 EU/mL, so

$$MVD = \frac{10 \ EU/mL \times 1}{0.05 \ EU/mL} = 200$$

## Product specific validation for testing FBS (fetal bovine serum) with the PyroMAT<sup>™</sup> system

### Assurance of the criteria for the endotoxin standard curve:

A standard curve using Reference Standard Endotoxin was performed to verify that the criteria for endotoxin standard curve were valid.

- The regression of response on log dose was statistically significant (p < 0.01)
- The regression of response on log dose did not deviate significantly from linearity (p > 0.05)

### Test for interfering factors and method validation for detection of non-endotoxin contaminants:

A dilution series from undiluted product up to the MVD was prepared and a test for interfering factors and method validation for non-endotoxin monocyte-activating contaminants (NEPs) according to EP was performed, comparing endotoxin-spiked sample dilutions as well as NEPs-spiked sample dilutions with the same unspiked sample dilutions.

Tests results showed non-endotoxin pyrogen (NEPs) control was detectable in all dilutions of the product from undiluted sample to MVD.

The undiluted product was found to show interference with the detection of the endotoxin spike, the reaction was strongly inhibited. From an 1:10 dilution to the MVD, the product reproducibly showed both detection of the NEP and a spike recovery within the range 50-200%, allowing to rule out interference with the test.

### Test for interference in the detection system:

The so found optimum dilution was forwarded to a test for interference in the detection system. A dilution series of IL-6 control was tested in absence and presence of the sample dilution. All IL-6 control dilutions showed less than 20% difference between dilution in absence and presence of the product, therefore no interference with the detection system could be found.

#### Conclusion

Assurance of criteria for the standard curve	Valid
Test for interfering factors	Valid
Detection of NEP contaminants	Valid
Interference in the detection system	Valid

All criteria of the product specific validation were fulfilled and the dilution 1:10 was chosen as the first valid dilution for the Method A.



### Quantitative method A according to EP

### **Planning Test Execution**

Once a valid dilution has been identified through the Product Specific Validation, Method A (quantitative test) can be performed to assess the pyrogenicity of samples from this product.

The lowest dilution / highest concentration of the sample, that was found free of interference in the PSV is to be chosen as initial dilution and 2-fold serial dilutions are tested.

For quantification, the middle section of the doseresponse curve of the endotoxin standard is to prefer, giving the most exact results, while quantification within the upper plateau of the standard curve can lead to inaccuracy due to reaching the endpoint of the reaction and therefore is not recommended.

For most exact results, we therefore recommend to analyze sample dilutions which do not exceed the measuring range of 0.05 to 0.4 EU/mL.

#### **Data interpretation**

The data analysis was performed with Gen5 software and the PyroMAT<sup>™</sup> Method A available on our website. Information related to the sample were completed directly on the software: sample name, CLC.

For interpretation, the layout for method A was modified with the appropriate dilution factors for this sample matrix.

latrix	Statistics	Graphs										-1-
Data: Plate Layout   Edit Matrix Read #: Show												
	1	2	3	4	5	6	7	8	9	10	11	12
А	BLK	BLK	BLK	BLK	SPL1:1 10	SPL1:1 10	SPL1:1 10	SPL1:1 10	SPL2:1 10	SPL2:1 10	SPL2:1 10	SPL2:1 10
в	STD1	STD1	STD1	STD1	SPL1_R	SPL1_R	SPL1_R	SPL1_R	SPL2_R	SPL2_R	SPL2_R	SPL2_R
D	0,0125	0,0125	0,0125	0,0125	10	10	10	10	10	10	10	10
С	STD2	STD2	STD2	STD2	SPL1:2	SPL1:2	SPL1:2	SPL1:2	SPL2:2	SPL2:2	SPL2:2	SPL2:2
Ŭ	0,025	0,025	0,025	0,025	20	20	20	20	20	20	20	20
D	STD3	STD3	STD3	STD3	SPL1_R	SPL1_R	SPL1_R	SPL1_R	SPL2_R	SPL2_R	SPL2_R	SPL2_R
	0,05	0,05	0,05	0,05	20	20	20	20	20	20	20	20
E	STD4	STD4	STD4	STD4	SPL1:3	SPL1:3	SPL1:3	SPL1:3	SPL2:3	SPL2:3	SPL2:3	SPL2:3
	0,1	0,1	0,1	0,1	40	40	40	40	40	40	40	40
F	STD5	STD5	STD5	STD5	SPL1_R	SPL1_R	SPL1_R	SPL1_R	SPL2_R	SPL2_R	SPL2_R	SPL2_R
<u> </u>	0,2	0,2	0,2	0,2	40	40	40	40	40	40	40	40
G	STD6	STD6	STD6	STD6	SPL1_NEF	SPL1_NEF	SPL1_NEF	SPL1_NEF	SPL2_NEF	SPL2_NEF	SPL2_NEF	SPL2_NE
	0,4	0,4	0,4	0,4	10	10	10	10	10	10	10	10
н	STD7	STD7	STD7	STD7	NEP	NEP	NEP	NEP				
	0,8	0,8	0,8	0,8	1	1	1	1				

### Testing of FBS (fetal bovine serum) with method A

The test setup was performed according to the user guide of the PyroMAT<sup>TM</sup> system.

An endotoxin standard curve was performed for the test.

Three dilutions of the product were tested according to method A described in EP: the dilutions 1:10, 1:20 and 1:40 of the sample were all tested with and without endotoxin spike. A quantification of the found endotoxin equivalents for all dilutions of the unspiked and spiked sample was done using the endotoxin standard curve. The endotoxin spike recovery for all three sample dilutions was calculated.


Heat killed *Staphylococcus aureus* (HKSA) was used as additional control for detection of non-endotoxin pyrogens within the sample, tested with the highest concentration of the product to be examined.

5

### After reading the plate, the data interpretation was performed with the software.

atrix	Statistics	Graphs										4
Data: Raw Data - Delta OD Summary (1/2) V 🕅 Edit Matrix Read #: 👘 Show												
	1	2	3	4	-5	6	7	8	9	10	11	12
А	BLK	BLK	BLK	BLK	SPL1:1	SPL1:1	SPL1:1	SPL1:1	SPL2:1	SPL2:1	SPL2:1	SPL2:1
~	0,013	0,013	0,016	0,016	0,017	0,017	0,019	0,018	0,000	0,000	0,000	0,000
в	STD1	STD1	STD1	STD1	SPL1_R	SPL1_R	SPL1_R	SPL1_R	SPL2_R	SPL2_R	SPL2_R	SPL2_R.
	0,013	0,015	0,015	0,017	0,242	0,259	0,311	0,288	0,000	0,000	0,000	0,000
с	STD2	STD2	STD2	STD2	SPL1:2	SPL1:2	SPL1:2	SPL1:2	SPL2:2	SPL2:2	SPL2:2	SPL2:2
Ŭ	0,025	0,018	0,017	0,020	0,017	0,016	0,015	0,015	0,000	0,000	0,000	0,000
D	STD3	STD3	STD3	STD3	SPL1_R	SPL1_R	SPL1_R	SPL1_R	SPL2_R	SPL2_R	SPL2_R	SPL2_R.
	0,057	0,042	0,034	0,024	0,245	0,255	0,230	0,298	0,000	0,000	0,000	0,000
F	STD4	STD4	STD4	STD4	SPL1:3	SPL1:3	SPL1:3	SPL1:3	SPL2:3	SPL2:3	SPL2:3	SPL2:3
	0,189	0,165	0,150	0,095	0,015	0,015	0,013	0,014	0,000	0,000	0,000	0,000
F	STD5	STD5	STD5	STD5	SPL1_R	SPL1_R	SPL1_R	SPL1_R	SPL2_R	SPL2_R	SPL2_R	SPL2_R.
1 - C	0,684	0,575	0,454	0,503	0,238	0,223	0,224	0,236	0,000	0,000	0,000	0,000
G	STD6	STD6	STD6	STD6	SPL1_NEF	SPL1_NEP	SPL1_NEP	SPL1_NEF	SPL2_NEF	SPL2_NEP	SPL2_NEP	SPL2_NE
Ŭ	3,216	2,735	2,815	1,750	1,809	1,895	1,698	1,946	0,000	0,000	0,000	0,000
н	STD7	STD7	STD7	STD7	NEP	NEP	NEP	NEP				
	4,136	4,308	4,314	4,355	1,894	1,615	1,504	1,691				
Н	STD7	STD7	STD7	STD7	NEP	NEP	NEP	NEP	0,000	0,000	0,000	

The standard curve was valid for all the criteria.

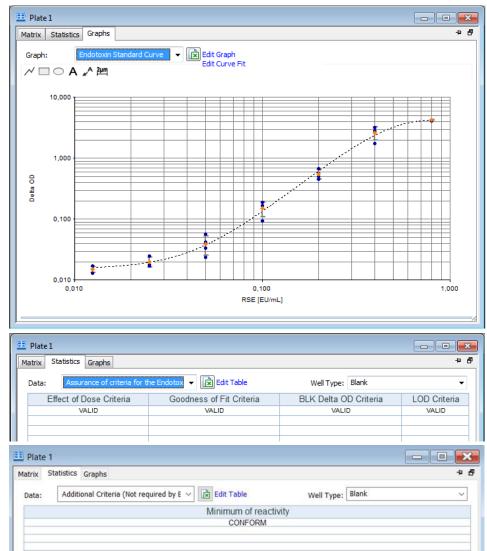

### The NEP-control confirmed detection of non-endotoxin pyrogens in the system and in the sample.

😐 Plate	1									×
Matrix	Statistics	Graphs								4 B
Data:	Dete	ction of No	on Endotoxin Py	rogen: 🔻	Edit Table	w	/ell Type:	Blank		•
			Dete	ection of I	Non Endotoxin	Pyrogens by the	e system	ı		
					VALI	D				
r										
🔠 Plate	1								_ 0	<b>x</b>
Matrix	1 Statistics	Graphs								× -
	Statistics		on Endotoxin Py	rogen: 🔻	Edit Table	W	/ell Type:	Blank		+
Matrix	Statistics					w n Pyrogens in Sa		Blank		
Matrix	Statistics					n Pyrogens in Sa		Blank		
Matrix	Statistics				f Non Endotoxir	n Pyrogens in Sa		Blank		
Matrix	Statistics				f Non Endotoxir	n Pyrogens in Sa		Blank		

The chosen sample dilutions were appropriate for the sample with all tested dilutions being within the MVD and spike recovery within 50-200%.

The sample itself showed a pyrogenic load of <0.5 EEU\*/mL which is below the CLC (10 EU/mL) and therefore being considered "not pyrogenic".

The data analysis showed validity of the test being fulfilled for all sample dilutions.

\* EEU: Endotoxin Equivalent Unit

Plate 1	L										_ 0
1atrix	Statistics	Graphs									÷.
Data: Data Analysis of Sample 1											-
Well I	D Name	Dil Factor	MVD	Delta OD	Mean	CV (%)	CLC [EU/mL]	[EEU/mL] x	DIL	Spike Rec. (%	6) CONCLU
SPL1	FBS	10	200,0	0,017	0,018	5,4	10,0	<0,500		69,3	PASS
		10		0,017							
		10		0,019							
		10		0,018							
		20	200,0	0,017	0,016	6,1	10,0	<1,000		67,3	PASS
		20		0,016							
		20		0,015							
		20		0,015							
		40	200,0	0,015	0,014	6,7	10,0	<2,000		64,1	PASS
		40		0,015							
		40		0,013							
		40		0,014							

### Results

The examined FBS showed inhibition of the reaction of the Monocytes in the undiluted sample that could be overcome by dilution within the authorized dilution range (not exceeding the MVD).

The data show that PyroMAT<sup>™</sup> system is suitable for detecting both endotoxin and non-endotoxin pyrogenic contaminations in preparations of FBS.

7

### References

[1] European Pharmacopoeia, chapter Monocyte activation test (2.6.30)

[2] Kirikae T, Tamura H, Hashizume M, Kirikae F, Uemura Y, Tanaka S, Yokochi T, Nakano M. "Endotoxin contamination in fetal bovine serum and its influence on tumor necrosis factor production by macrophagelike cells J774.1 cultured in the presence of the serum" International Journal of Immunopharmacology Volume 19, Issue 5, May 1997, Pages 255-262

[3] Fishel S, Jackson P, Webster J, Faratian B. "Endotoxins in culture medium for human *in vitro* fertilization" Fertility and Sterility Volume 49, Issue 1, January 1988, Pages 108-111

[4] Ryan J. "Endotoxin and cell culture" Technical bulletin, www.corning.com/lifesciences



MilliporeSigma 400 Summit Drive Burlington, MA 01803

#### To place an order or receive technical assistance

In the U.S. and Canada, call toll-free 1-800-645-5476

For other countries across Europe and the world, please visit: **EMDMillipore.com/offices** 

For Technical Service, please visit: EMDMillipore.com/techservice

#### EMDMillipore.com/pyromat

MilliporeSigma, the vibrant M, Millipore and PyroMAT are trademarks of Merck KGaA, Darmstadt, Germany or its afiliates. All other trademarks are the property of their respective owners. Detailed information on trademarks is available via publicly accessible resources. © 2018 Merck KGaA, Darmstadt, Germany and/or its affiliates. All Rights Reserved.

MS\_AN1804EN Ver. 2 2018-10080 04/2018