Supelco

Analytical Products

Reflectoquant®

RQflex® 20

Quick Guide

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Guide rapide

Guía Rápida

Guia Rápido

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EN

Quick Guide

For the Operating Manual please visit:
www.sigmaaldrich.com/rqflex

EN

Table of contents

1	Intended use	5
2	Getting started	5
2.1	Package contents	5
2.2	Inserting the batteries	6
2.3	Operation buttons and Display	7
2.3.1	Operation buttons	7
2.3.2	Display	8
2.3.3	Menu items	9
2.4	Settings	9
2.5	Initial calibration	10
3	Measurement	11
3.1	Add a new method	11
3.2	Measuring with test kits	13
3.3	Notes on the measurement	15
3.4	Application notes	16
4	Maintenance of the instrument....	16
4.1	Handling	16
4.2	Cleaning the strip adapter	17
5	Troubleshooting	17
6	Technical data	18
7	Accessories	18
8	Service / Warranty	19

1 Intended use

The Reflectoquant® RQflex® 20 reflectometer is intended to be used exclusively according to the operating manual. The reflectometer was developed for performing analyses of chemical parameters in water, food & beverage samples, and environmental samples using Reflectoquant® test strips in the laboratory, on the production line or in the field. Any other use is considered to be unauthorized.

EN

2 Getting started

2.1 Package contents


The standard contents of the Reflectoquant® RQflex® 20 reflectometer package comprise the following items:

- 1 Reflectometer including strip adapter
- 1 Recalibration set
- 1 Quick Guide
- 4 1.5-V batteries (AAA)
- 1 Safety Instructions

2.2 Inserting the batteries

Before operating the system for the first time, the batteries included in the package must be installed.

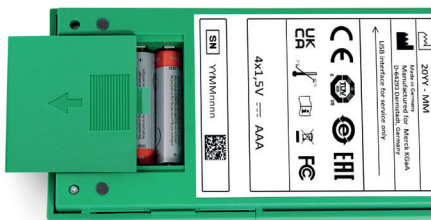Dispose of used batteries in accordance with the local regulations.

Insert the batteries into the compartment, heeding the + and – pole indicators.

Recommendation

Do not use rechargeable batteries!

Saving data - Important notes


Before changing the batteries, please shut down the instrument first. Otherwise the date and time settings are lost. During battery change, data in the RQflex[®] 20 reflectometer is saved for 60 seconds. If the change time exceeds 60 seconds, only date and time are lost, all stored data and settings are still stored.

Recommendation

Replacement batteries (AAA non-rechargeable alkaline batteries, 1.5 V) must be on hand for immediate insertion.


2.3 Operation buttons and Display

2.3.1 Operation buttons

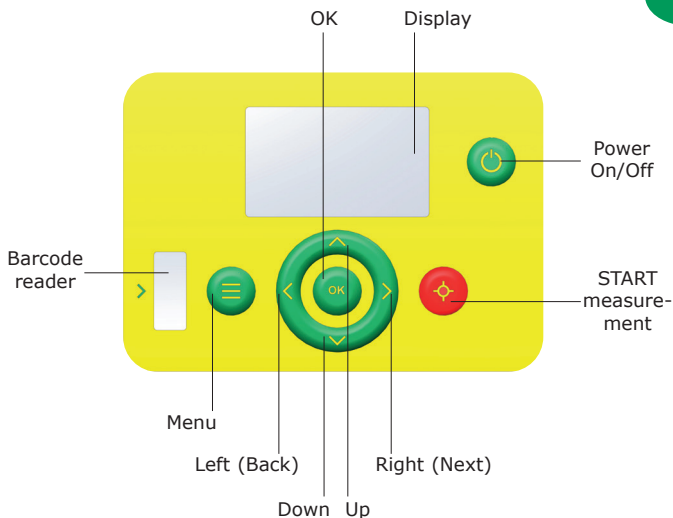
2.3.2 Display

General Structure

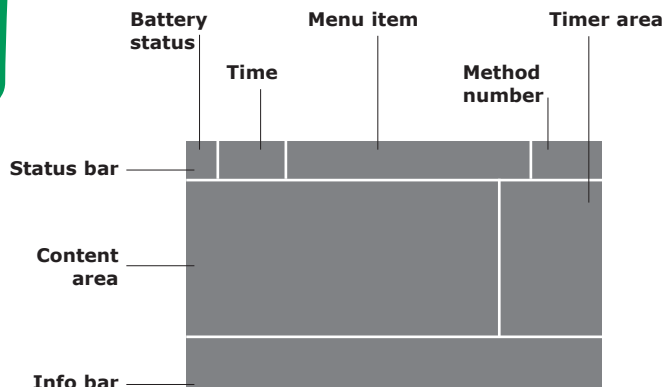Example

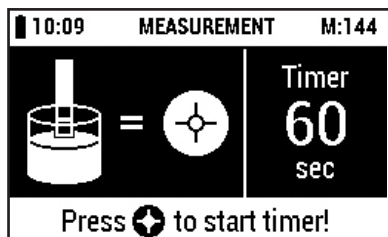2.3.3 Menu items

The following items are displayed in the main menu:

Method List	list of all stored methods
Result List	list of all stored results
Settings	list of instrument settings
Quality Assurance	list of all analytical quality assurance options
Information	list of instrument information

The logo consists of the letters 'EN' in a bold, yellow, sans-serif font, positioned within a green, rounded rectangular background.

2.4 Settings


The following items are displayed in the settings menu:

Languages	for setting the desired language (EN - DE - FR - ES - PT - JP - CN)
Date Format	for setting the desired date format (yyyy-mm-dd - dd.mm.yyyy - mm/dd/yyyy)
Time format	for setting the desired time format (24 h - am/pm)
Date	for setting the current date
Time	for setting the current time
Auto Power Off Time	for setting the desired power off time (0 (= Off) - 15 - 30 - 45 - 60 min)
Audible Signal	for setting the audible alert (On / Off)
Acoustic. Countdown Warning	for setting the desired duration of the warning signal for the countdown (0 (= Off) - 5 - 10 - 15 - 20 - 25 - 30 sec)
Memory Capacity Warning	for setting the memory capacity warning (On / Off)


2.5 Initial calibration

The instrument needs to be calibrated before the first measurement is performed. This is required to set the correct baseline for the reflectometric optics. For this purpose, a Reflectoquant® Recalibration Set is delivered with the instrument. Please use the barcode and the plastic test strip for calibration (for details of the Recalibration Set see Operating Manual). Keep the instrument at ambient temperature for at least 30 minutes.

Select [**Quality Assurance**] in the [**MAIN MENU**].

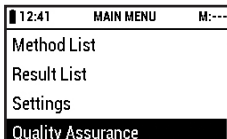Select [**Calibration**].

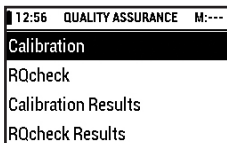Confirm your selection by pressing [OK] or [Right (Next)].

Follow the instruction shown on the display.

Correct insertion of the test strip:

Insert the calibration strip of the recalibration set into the strip adapter with the beveled edge outwards and upwards.

A **calibration** has to be repeated

- after changing/cleaning of the strip adapter and/or the internal standard
- in the case of issues during the boot up (self-check test failed)
- in the case of error messages E07 or E10 (after cleaning the optics)
- in the case of RQcheck failure (error messages E09 and E11)
- in the case of suspect measurement results
- after severe mechanical distress (e.g. after dropping the instrument).

3 Measurement

There are two scenarios on how to start a measurement:

1. The desired method is not yet saved in the instrument.
Continue with chapter 3.1.
2. The method is already saved in the instrument.
Continue with chapter 3.2.

3.1 Add a new method


There are two scenarios for adding a method:

1. Guided way
2. Shortcut


Guided way

EN

Select [**Method List**] in the [MAIN MENU].

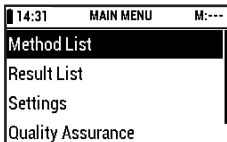Select [**Add new Method**].

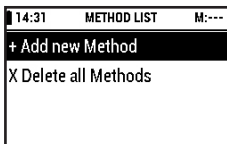
Confirm your selection by pressing [OK] or [Right (Next)].

Follow the instruction shown on the display.

The display shows the new added method in the first row.

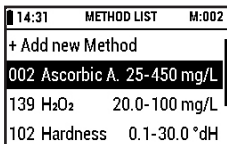Remove the barcode strip from the barcode reader and store it in the supplied packaging material of the Reflectoquant® test (not in the test-strip tube).

Press the [OK] button for method details.

14:31	METHOD DETAIL	M:002
Article Number 1.16981.0001		
Method	Ascorbic A. 002	
Range	25-450 mg/L	
X Delete Method		

Shortcut

A method can be saved and started directly by inserting the barcode strip of the Reflectoquant® test.

3.2 Measuring with test kits

There are three ways to measure with test kits:

1. Continue with last used method
2. Choose from method list
3. Shortcut

Continue with last used method

The method number shows the last used method, measurement is started immediately after pressing the red [START] button.

12:38	MAIN MENU	M:002
Method List		
Result List		
Settings		
Quality Assurance		

Choose from method list

If another stored method is needed, select [**Method List**] in the [**MAIN MENU**].

12:38	MAIN MENU	M:002
Method List		
Result List		
Settings		
Quality Assurance		

Select the desired method.

12:39	METHOD LIST	M:002
+ Add new Method		
029 K	0.25-1.20 g/L	
002 Ascorbic A.	25-450 mg/L	
102 Hardness	0.1-30.0 °dH	

Confirm by pressing the [START] button

and follow the instruction in the package insert of the corresponding Reflectoquant® test, e. g. "Preparation", "Procedure".

Shortcut

A method can be started directly from any menu item by inserting the barcode strip of the Reflectoquant® test into the instrument.

Follow the instruction in the package insert of the corresponding Reflectoquant® test, e. g. "Preparation", "Procedure".

Note



For a detailed measurement workflow description, see Operating Manual.

3.3 Notes on the measurement

Code number

The first number (method number) corresponds to the code number of the barcode strip, or respectively the first three digits of the batch number of the Reflectoquant® test strip package:

22:29	METHOD LIST	M:002
+ Add new Method		
002	Ascorbic A. 25-450 mg/L	
139	H ₂ O ₂	20.0-100 mg/L
102	Hardness	0.1-30.0 °dH

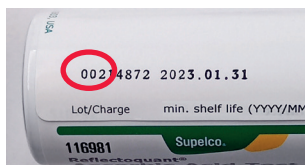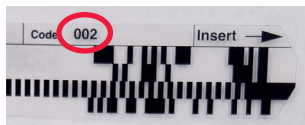Skip the timer

It is possible to cancel the timer countdown, to perform measurements directly without the countdown function of the instrument. Still, make sure that the reaction time of the test method is met.

For details, see Operating Manual ("Skip the timer").

Serial measurements

It is possible to interrupt the timer to perform serial measurements.

For details, see Operating Manual ("Serial measurements").

3.4 Application notes

Various detailed application notes are available to support your analysis, especially the sample preparation step of your workflow:

www.sigmaaldrich.com/wfa-applications

4 Maintenance of the instrument

4.1 Handling

Please treat this instrument with the same care as you do with all other electronic devices. Ensure that no liquid enters the instrument, since this may result in damage that is not covered by the warranty.

Do not expose the instrument to excessively humid conditions or to excessive heat or cold (see Technical Data).

If the external casing of the instrument needs cleaning, please wipe the casing and the display carefully with a moist cloth.

Clean the optics and the strip adaptor only with water and a mild detergent or ethanol, 70% (v/v) on a tissue or cloth after dismantling the strip adaptor (see 4.2). Do not rinse the instrument or parts of it.


4.2 Cleaning the strip adapter

The strip adapter should be thoroughly cleaned (dismantling the strip adapter into its individual components) at the end of each working day.

The strip adapter must be cleaned immediately if an error message (see Operating Manual, "Troubleshooting") is displayed.

Clean the strip adapter only with distilled water and a mild detergent or ethanol, 70% (v/v) on a tissue or cloth after dismantling the strip adapter. Do not rinse the instrument or parts of it.

EN

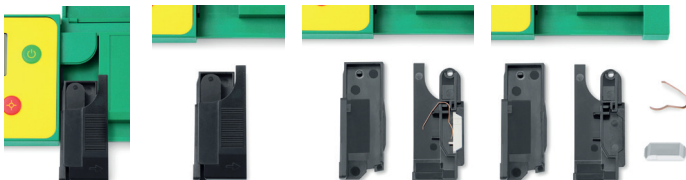5 Troubleshooting

The most commonly encountered problems are caused by

- incorrect insertion of the test strip
- incorrect reaction time
- incorrect use of the tests (e. g. strips not properly stored, pH range not correctly adjusted).

For more details, see the Operating Manual ("Troubleshooting"). The tables there show explanations of the error messages and tips on how to avoid incorrect measurements.

6 Technical data

Dimensions:	184 x 79 x 30 mm
Weight:	253 g (including batteries)
Memory:	50 test methods, 200 measurement results
Interface:	yes (for technical service only)
Light source:	4 LEDs, green/red, double optics
Power source:	4 x 1.5-V batteries (AAA)
Display:	Reflective LCD Graphic Module (256 x 160 dot)
System diagnosis:	yes
Measurement range:	4 - 90% rel. remission
Reflection area:	4 x 6 mm
Resolution:	0.1% rel. remission
Operating temperature:	5 - 40 °C for ideal measurements
Operating humidity:	below 80% for ideal measurements

7 Accessories

Reflectoquant® Strip adapter Replacement part for RQflex® 20	1.17267.0001
Reflectoquant® Recalibration Set Replacement part for RQflex®	1.16954.0001
Reflectoquant® RQcheck set for RQflex® 20 Accessory for RQflex® 20	1.17247.0001

8 Service / Warranty

Our instruments are 100% quality-controlled, which means that each individual instrument has been tested before leaving the factory. In addition to the described cleaning and calibration of RQflex® 20, the device is maintenance-free. The detailed instruction should allow for problem-free operation.

In the event that you encounter an unsolvable problem, please contact the technical support specialist.

www.sigma-aldrich.com/customer-service.html

EN

Warranty

The manufacturer grants for this RQflex® 20 reflectometer a total replacement warranty of 12 months from the date of purchase.

All warranty work is provided by our customer service unit. In the event of a proven production or material defect, we will provide you with a new or as-new replacement RQflex® 20 reflectometer free of charge.

Improper handling results in loss of warranty. No water or any other liquid must be allowed to enter the instrument. In such cases, the warranty becomes invalid.

All warranty rights become invalid in the event that the purchaser or an unauthorized third person repairs or opens the instrument, or if changes are made to the warranty slip. Only the customer service unit is authorized to carry out repairs on components.

If the instrument is sent in for repair under the terms of the warranty, a copy of the invoice or of another proof of purchase must be enclosed in all cases.

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The logo consists of the letters 'PT' in a bold, yellow, sans-serif font. The letters are positioned on a green background that has a rounded, irregular shape, resembling a stylized leaf or a piece of paper. The 'P' and 'T' are closely spaced.

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The logo consists of the letters 'PT' in a bold, yellow, sans-serif font, positioned on a green background that is part of a larger graphic element on the right side of the page.

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
We provide information and advice to our customers on application technologies and regulatory matters to the best of our knowledge and ability, but without obligation or liability. Existing laws and regulations are to be observed in all cases by our customers. This also applies in respect to any rights of third parties. Our information and advice do not relieve our customers of their own responsibility for checking the suitability of our products for the envisaged purpose.

The life science business of Merck operates as MilliporeSigma in the U.S. and Canada.

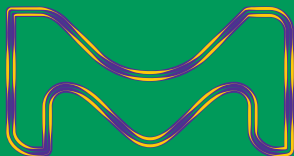Merck KGaA, 64271 Darmstadt, Germany
Tel. +49(0)6151 72-2440

www.sigmaaldrich.com/rqflex

Merck, the vibrant M, Supelco, Reflectoquant, and RQflex are trademarks of Merck KGaA, Darmstadt, Germany or its affiliates. All other trademarks are the property of their respective owners. Detailed information on trademarks is available via publicly accessible resources.

© 2024 Merck KGaA, Darmstadt, Germany and/or its affiliates. All Rights Reserved.

MERCK